

Softband 5

Variety of sizes and colors

The perfect options for just the right fit

Two adjustable sizes

- XS/S: 15-20 inches (380-500mm) head circumference
- M/L: 20-24 inches (500-600mm) head circumference

Unilateral and bilateral versions




Unilateral



Bilateral

In 8 different colors

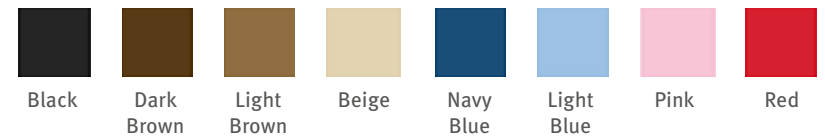

Providing impressive sound transmission

Newly designed connector plate and connector pad allow for greater sound transmission



New connector pad



- Unilateral softband contains 3 connector pads
- Bilateral softband contains 6 connector pads
- Replacement packs contain 20 pads



Softband 5 without pad	Softband 5 with pad
A top-down view of the softband 5 without the connector pad, showing the black fabric and the metal connector plate.	A top-down view of the softband 5 with the connector pad, showing the black fabric, the metal connector plate, and the black cushioned pad.
A side view of the softband 5 without the connector pad, showing the black fabric and the metal connector plate.	A side view of the softband 5 with the connector pad, showing the black fabric, the metal connector plate, and the black cushioned pad.
A heatmap showing the force distribution of the connector plate when fitted on a patient's head without the connector pad. The force is concentrated in a small area, with a peak of 121N.	A heatmap showing the force distribution of the connector plate when fitted on a patient's head with the connector pad. The force is more evenly distributed, with a peak of 23N.
<p>These are measurements of force distribution of the connector plate when fitted on a patient's head with and without the connector pad, measured in Newton.</p> <p>The connector pad helps provide a more targeted and efficient transfer of the sound signal.</p>	